

NO.11

MZ522

A magical fact!
An amazing effect!



产品代言人:


国际巨星 Aimor 女士

2ND GENERATION INTELLIGENT SKIN REJUVENATOR

★ ★ ★ *Product Instructions* ★ ★ ★

Contents

Foreword	1
Operating Principle	2
Suitable Situations	4
Structure and Performance	4
Directions	4
Attentions	5
Appendix list	5

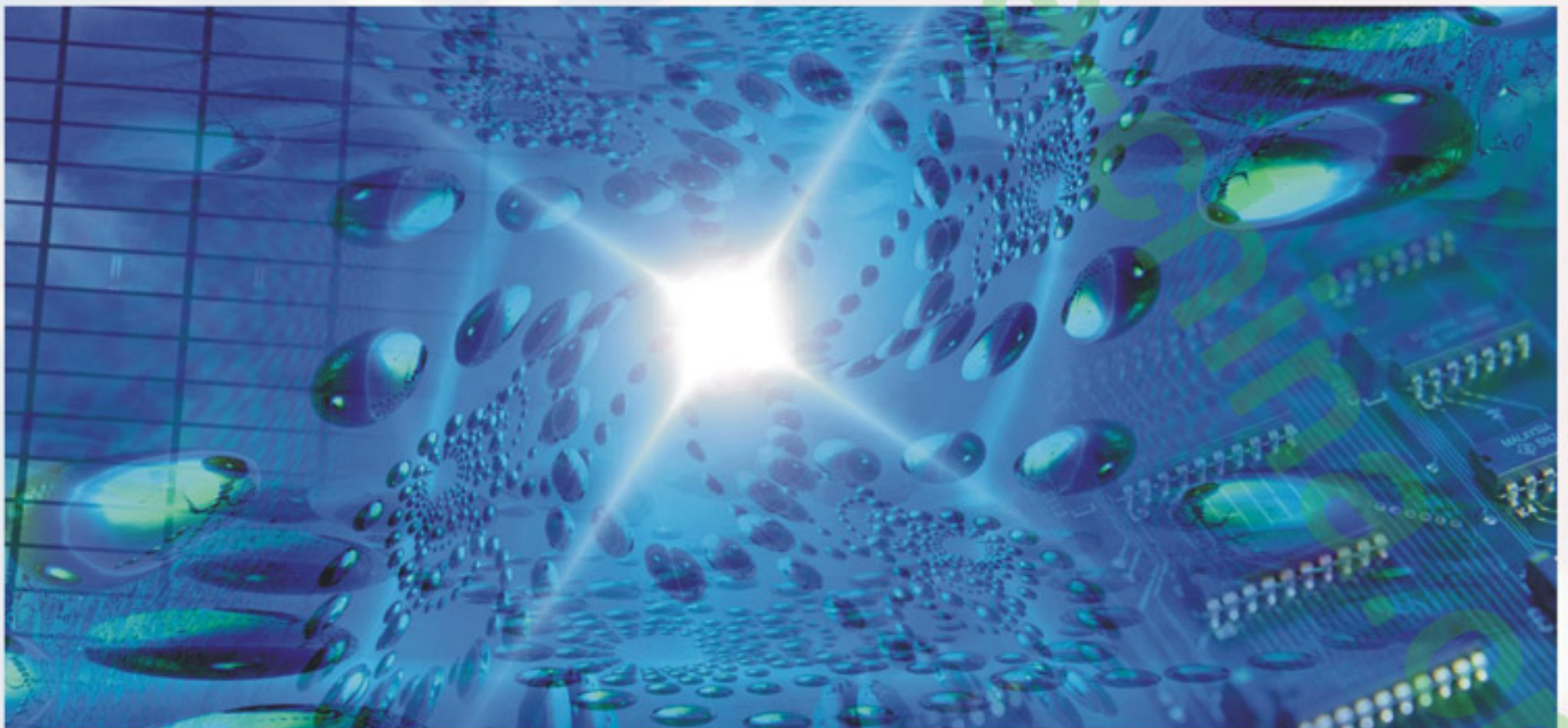

Foreword

Once, one of my friends said to me, “There used to be a youthful figure for me, but I didn't cherish it. And I was so regretful when it was gone. It is the largest pain in my life. If the God would give me another chance, I would say: I need youth. If need to add an deadline, I hope it is forever”. Imitating a classic dialogue from “A Chinese Odyssey” by Stephen Chow, my friend was quite humorous. But her words indicated some sadness and helplessness.

Youth is a priceless treasure. Life is only once for everybody. Within one life, the most colorful years should be the short youthful period. River flows and time flies. The marks of age emerge on our face cruelly. You may forget the joy when your height finally increased by one inch, but you would never forget the shock wrinkles on your face at the first time. You may forget the satisfaction to a fresh summer and dress the new skirts, but how can you forget your facial spots resulted from life and work pressure?

Life is loaded with too many dreams. But everybody has the same dream. That is: I want to be youthful again!

Now, we herein make an announcement to the public: the “**2nd Generation Intelligent Skin Rejuvenator**” is now available! To be youthful again is no longer a dream!

For youthful girls, aren't you afraid of meeting your boy friends with acnes on your face? For successful ladies, can't you let the marks of age mount onto your angles of eyes? For handsome gentlemen, don't you mind your dull face that overshadows your shiny heart?

NO! NO! The answer is always “NO”! So lucky people, please pay attention to the super youthful show performed by the “**2nd Generation Intelligent Skin Rejuvenator**” in China. If you're bothered by wrinkles, acquired color spots and color dots, or acnes, pimple, pouches, dark circles, dullness or roughness, the Rejuvenator will be helpful! This “beauty angel” from America is equipped with a wonderful operating system and a high-technology, intelligent and multifunctional repair network, to help you restore a charming and youthful appearance.

Besides remarkable benefits, this Rejuvenator is characterized by simple, convenient and time-saving operation. Just put the mask on face gently, tie the motor processor to your waist, and press the starting button. Within 16-minute beautification, you can do whatever you want, like reading a book, watching TV, chatting with friends, or even doing housework... You may remember your holidays when you spent at least one hour sitting or lying silently in a beauty parlor. Now you know beauty care is so easy; how can you resist the charm of this perfect skin rejuvenator?

Operating Principle

The “2nd Generation Intelligent Skin Rejuvenator” is armed with Cell Aging Control Technology, an up-to-date technical achievement in the 21st century. It solves the problem of skin aging fundamentally. Its rejuvenating repair system is composed of an RF facial mask, a main processor, Rejuvenator collaboration Series (1,2,3,4,5,6,7) seven parts and IPL rejuvenation Instrument (mainly used laser beyond the advanced technology of optical pulse on the medical). It contains a variety of different wavelengths in the spectrum of skin from different depths within the regional temperature changes, optical effects, particle vibration damage of the increase in bad cells, amino acids, proteins and other life material to enhance the activity of the skin to repair skin metabolic effects. Composite optical pulse technology is a non-injury treatment technology breakthroughs can reach the depths of the skin treatment regional. the result is very security and reliable). Among others, its motor processor equivalent to a Pentium 7 processor in computer is playing a critical role. Intelligent chip of the processor is acting like a brain. For different skin conditions, the processor can formulate combined energy and deliver it to the RF facial mask by means of multi-guided wave. Together with 24 gold-plated inductive points on the RF facial mask, it jointly forms an intelligent repair network. Like a “beauty doctor”, this network can diagnose and position damaged skin cells automatically. Aiming at the damaged cells, the network can transmit corresponding combined energy to repair skin intelligently. Different users may have individualized feeling during repair, which reveals circumspection of the skin rejuvenator.

This Rejuvenator is developed on the basis of physics, electrochemistry, dialysis, iontophoresis and relevant theories. Cell is the basic unit of human body. Cell membrane controls the exchange of substances inside and outside the cell. Utilizing absorptive and repulsive forces of electric field, the cell membrane can absorb essential nutrients into the cell. Based on this theory, this Rejuvenator adopts up-to-date electronic and bioscience technology from America to realize perfect skin care.

What magical changes does the “2nd Generation Intelligent Skin Rejuvenator” bring to customers?

The results of follow-up survey on 100,000 users clearly show that:
After 16-minute application, their skin become smooth and youthful.
After 3-day application, their skin turn up exquisite, white and shiny.
After 1-week application, their facial fine lines and dark spots start to decrease; their pouches and dark circles begin to disappear.
After 4-week application, they're feeling 1-2 years younger.
After over 1-year application, they look 5-10 years younger than their ages.



Then, how come the amazing rejuvenating effects?

Uptodate researches in the world proved that: the deliverability of chondriosome in cell decides cell composition and division. Once the deliverability is weakened, the speed of cell division and renewal will be slowed down, which accordingly leads to reduction of new cells and increasing of aged cells, as well as appearance of wrinkles, looseness and pouches.

This Rejuvenator takes the lead to adopt biological waveguide intellectual technology, which settles the problem of insufficient cellular energy thoroughly. It helps eliminate wrinkles, pouches, looseness and sagginess, restoring youthful skin rapidly. The Time magazine evaluated it as the 4th top technology that affects American's life, an innovative achievement of modern beautification. It is a breakthrough of traditional skincare methods that only focus on epidermal absorption, cellular nourishing, skin regulation, etc.

The Rejuvenator emits activating biological signals that aim at removing wrinkles, looseness, sagginess and pouches resulted from skin aging. This is a breakthrough of traditional skincare methods that simply emphasize supplement of moisture and protein. The signals can act on deep cells directly to enhance the deliverability of chondriosome. Consequently, they accelerate cell division and renewal, promote the composition of collagen and elastic fiber, and maintain normal support of netted texture, eliminating skin aging fundamentally. The Rejuvenator emits regenerating biological signals that work to promote activity of enzymes in cell. They speed up intercellular biochemical reactions (DNA duplication and transcription) and promote cell division. Accordingly, they eliminate wrinkles, looseness, sagginess and pouches, leaving skin natural and youthful. Hydrolytic biological signals emitted by the Rejuvenator can act on withered, curled and abnormal cells directly, break them up rapidly, and discharge them by means of blood circulation and lymph circulation. Moreover, the Rejuvenator is beneficial to regulating excessive sebum, removing acnes and pimples thoroughly. Leaves skin youthfully smooth and elastic.

Circulating biological signals emitted by the Rejuvenator can act on deep cells of skin. Help promote metabolism and blood circulation, speed up nutrients absorption, and discharge metabolite waste rapidly. Keep skin active and free from dullness, roughness, dark circles and color spots. Leave skin white and natural, smooth and youthful.

Effective rates of the Rejuvenator:

- 100% in minimizing pores and refining skin
- 100% in eliminating acnes and pimples
- 88.9% in controlling wrinkles
- 90.9% in reducing and removing acquired color spots
- 87.5% in relieving or removing pouches and dark circles



Suitable Situations and Contraindications


1. Suitable Situations

Suitable for 25-60-year-old people in pursuit of beauty.

2. Prohibited for the following people:

- 1) People carrying pulse regulator or other electromagnetic devices;
- 2) People with serious heart disease;
- 3) Women during pregnancy and lactation period;
- 4) People who used to receive operation of heart transplantation;
- 5) Children under 6 years old.

Structure and Performance

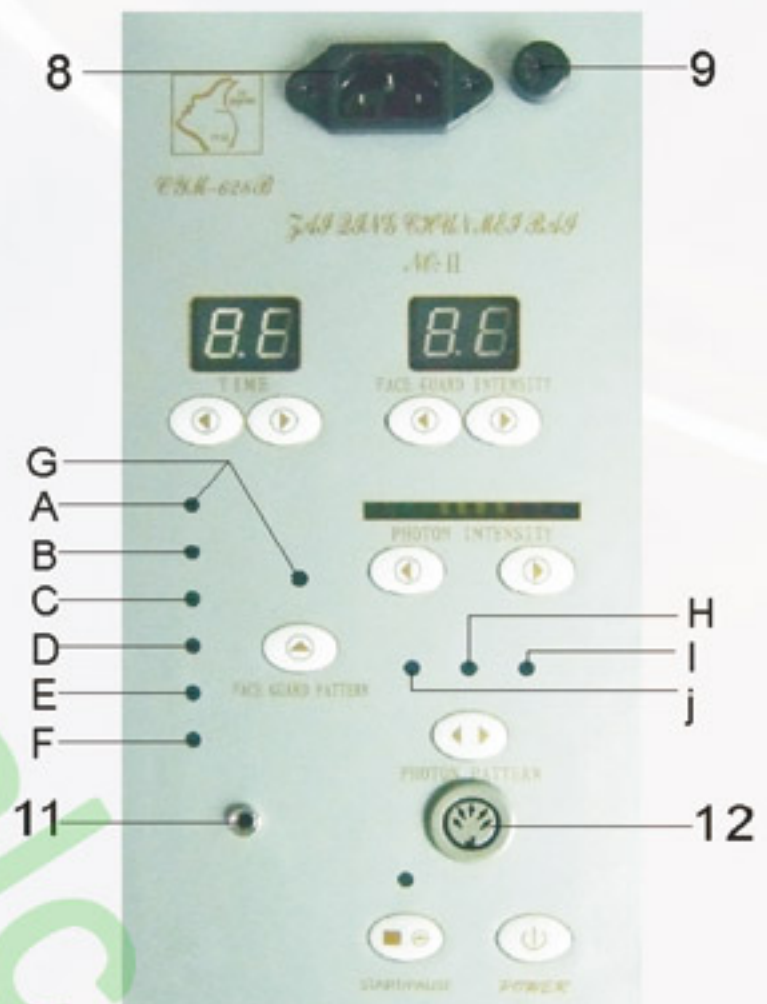

Directions

- A. Please check the accessories in line with the packing list, and read the Description before applications.
- B. Connect controller and photon facial mask of the rejuvenator, and fill in the specified battery.
- C. Apply the rejuvenating energy essence evenly to face according to skin conditions.
- D. Fix photonic band with a bandage mask on the face
- E. 1. Is the time choice key of photon mask, the time display show about 16 minutes when you turn on the machine, and it can be adjusted. (Generally used for 10 to 30 minutes)
2. Is the start button and select button of the photon masks, is divided into strength (0-30 files) can be different in accordance with our stage of the skin, so we can choose the strength
3. Is the 7 kind of working mode of the photon mask (a,b,c,d,e,f,g), we can choose different mode according to different skin:

- A. Intermittent waves----proud skin, diluted melanoma. ningbo delicate skin, a shinky.
- B. Low-frequency wave ----reduce the fat accumulation face and neck ,to enhance facial lines.
- C. Small sawtooth ---- reduce facial wrinkles: tight skin becomes soft, smooth, flexible.
- D.frequency wave ----wipe out wrinkles, facial skin acupoint massage.
- E.mutation wave ----reduce fat accumulation face and neck to eliminate double chin, facial lines upgrade.
- F.Middle frequency wave ----eliminate facial inflammatory acne: regulate sebum secretion, smooth the skin ,reduce the greasiness.
- G.Is for the whole cycle of computer memory automatically. (six working modes cycle automatically)
- 4.the intensity key of the photo, use to start IPL rejuvenation and adjust the photon intensity control button.
- 5.IPL rejuvenation instrument (H,I,J) has three different modes ,so we can select the key by different skin:h, 30 percent use optical frequency band and the promotion of collagen protein and activation of cell surface repair.
- I.And 50% for the optical frequency band contraction pores, skin tight it is bactericidal anti-inflammatory, healing sore face elimination and Downplay acne scars, effective prevention and Anchuang Jingchundou Growth
- J. 80% optical frequency band is able to dilute stain various age grouper face increased through the lymphatic and blood detoxification Central. (to use the IPL rejuvenation at the end of photond mask will get better results)
- 6.start Pause button. 7.On/off button . 8.power outlets. 9fuse blocks.
- 10. photon masks connectivity. 11. photon connectivity.

Precautions


1. After low voltage conversion, voltage of the device is safe. But soaking of controller is not allowed.
2. For fresh users, please start from the high gear and apply once every 3 days, 10 minutes once. After 5 applications, apply once a week.
3. After application, do not let water flow into the control panel.
4. Every time after application, please keep the rejuvenator in a dry place.
5. During the course of application, please drink more water as required. In case of other treatment, do not discontinue use at once.

List of Accessories

- Photon mask-----1pcs
- Photon -----1pcs
- Fuse-----1pcs
- Instruction-----1pcs
- Power line-----1pcs
- Photon mask connect line-----1pcs



***Fashionable in America,
Available in China now***

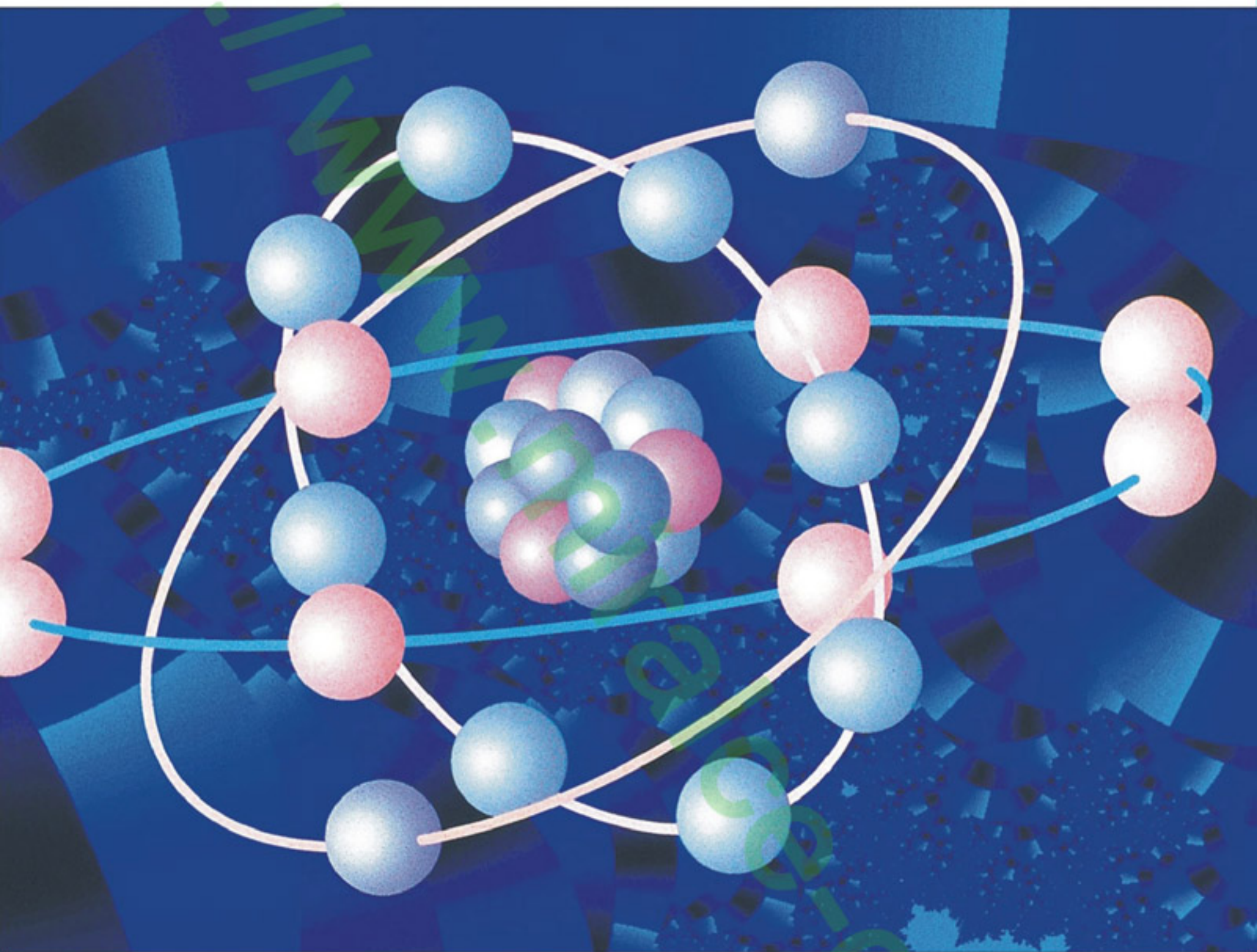

A secret to keep women's youth!

MIRACLE CO.,LTD

Address: Room B601 Hengsheng Building B , No. 205 Huanshi Zhong Road

Guangzhou China /POST NO.:510010

TEL:86-20-86673107 FAX:86-20-86672750

Email: admin@miracle-china.com

Skype: miracle-china

MSN:china.miracle@hotmail.com QQ:286372662

<http://www.miracle-china.cn>